



VERTUS

SEATING SYSTEM



For over 50 years, RMS (Rehabilitation Manufacturing Services) has been designing and manufacturing postural support and mobility equipment for children and adults with complex needs.

Trusted by NHS Wheelchair Services, therapists and healthcare professionals in schools and CDCs across the UK, we have built a reputation for dependable products, clinical expertise and long-term support.



A TRUSTED PARTNER FOR SCHOOLS AND CDCS

As a British manufacturer, we maintain full control over design and production, ensuring consistent quality, reliable supply and the flexibility to meet individual requirements. Our range of seating systems, mobility products and postural support accessories is developed to improve comfort, participation and functional outcomes across educational, clinical and home environments.

Through our nationwide team of Product Specialists, we work closely with schools, therapists and families to assess needs, recommend appropriate solutions and provide ongoing support. From initial assessment through to delivery and follow-up, our aim is to make the process as straightforward and effective as possible for both staff and users.

WHY RMS

- **50 years** of specialist manufacturing experience
- **Trusted by NHS and Schools** across the UK
- **British-made**, quality-controlled production
- **Nationwide** assessment and product support
- **ISO 13485 certified** for quality and safety

VRTUS



Made-to-Measure Support with Exceptional Adjustability

The Virtus Seating System is a highly adjustable, modular seating solution designed to support children and young adults with a wide range of postural needs within educational environments. Available in four sizes and individually configured to each user, Virtus combines adaptable support, growth adjustment and provides long-term seating that evolves alongside changing needs.

Designed around both clinical outcomes and everyday school life, the Virtus can be configured with either a fixed or dynamic backrest and paired with a manual or powered base, creating a flexible seating solution that supports learning, interaction and participation throughout the school day.

Key Features

- **Four sizes** supporting users from early years through to adulthood
- **Highly adjustable** design to accommodate growth and changing postural needs
- **Dynamic seating** technology to support movement and tone management
- Full range of **accessories**

Modular **Postural Support**

The Virtus Seating System is supplied with a comprehensive range of postural support accessories, enabling highly individual positioning and clinical configuration.

To accommodate the wide range of available options and customisations, Virtus seating systems are specified using our intuitive and interactive prescription form.

Trunk & Thoracic Support


- Fixed and swing-away laterals
- Shoulder and tricep supports
- Harness systems


Pelvic Positioning

- 2-point and 4-point belts
- Stabiliser belts
- Pelvic pads

Cushions & Covers

- Flat and contoured cushions
- Multiple foam options
- Washable and wipe-clean cover options

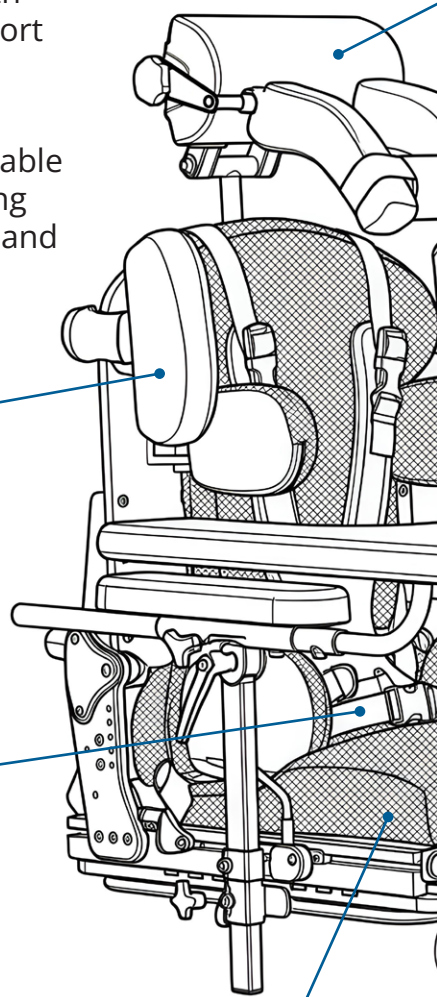

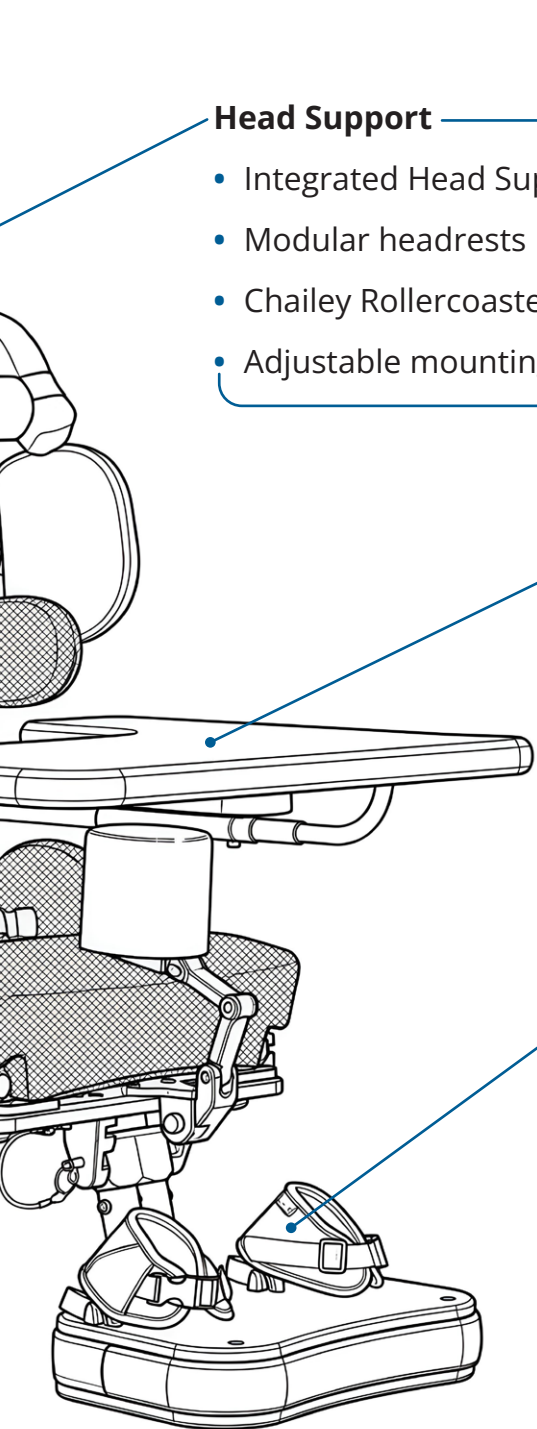

Head Support

- Integrated Head Supports
- Modular headrests
- Chailey Rollercoaster
- Adjustable mounting

Arm Support & Trays

- Adjustable arm supports
- Padded and moulded tray options
- Functional support for everyday activities


Leg & Foot Support

- Pommels
- Knee abductors
- Dynamic footrests
- Anklesures and foot control options



Fixed Backrest

The fixed backrest provides stable and effective postural support for users who do not present with significant involuntary extension or high tone.

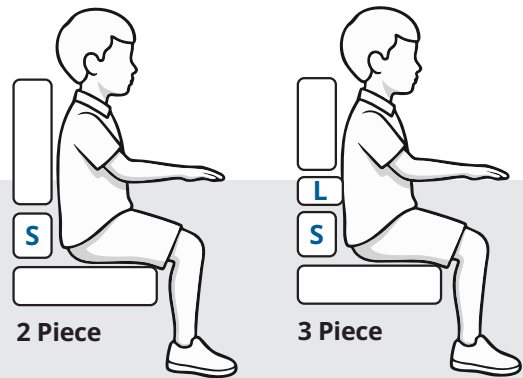

Available in both **two-piece** and **three-piece** configurations, each backrest is manufactured to the user's individual measurements. The overall backrest size, together with the dimensions of the sacral and lumbar support sections, can be tailored to provide highly personalised postural support.

	Size 1	Size 2	Size 3	Size 4
WHEELCHAIR SIZE	36	40	45	50
SEAT DEPTH	19 - 37 cm	26 - 45 cm	34 - 52 cm	34 - 52 cm
PELVIC WIDTH	21 - 32 cm	25 - 36 cm	30 - 41 cm	35 - 46 cm
THORACIC WIDTH	19 - 36 cm	25 - 40 cm	31 - 46 cm	35 - 51 cm
SACRAL PAD HEIGHT	10 cm	12.5 cm	15 cm	15 cm

Dynamic Support

For users who experience involuntary extension, fluctuating tone or forceful movements, the dynamic backrest allows controlled movement whilst maintaining postural support. By absorbing and managing these forces, the system can help improve comfort, seating tolerance and positioning throughout the day.

The dynamic mechanism is available with three resistance levels, allowing the response of the system to be matched to the user's individual needs. Low, medium and high resistance options provide flexibility across a wide range of postural presentations.



Long-Term Adjustability

The Virtus Seating System is designed to adapt alongside the user over time. Overall backrest height is individually configured using the user's shoulder height measurement with an additional 7.5cm growth allowance incorporated into the design.

Combined with extensive adjustment throughout the seating system, this approach helps maximise product lifespan whilst accommodating changing postural requirements.

VRTUS

⤴ HIGH - LOW ⤵



Max. Client Weight:

HI/LO BASE 90kg

MOBILITY BASE 135kg

Powered Positioning with Everyday Practicality


The Virtus Powered High-Low Base is a robust and highly adaptable indoor seating base designed to support positioning, interaction and participation throughout the school day. Combining powered height adjustment and tilt functionality with a durable heavy-duty frame, the system allows seating positions to be adjusted easily for learning, therapy, transfers and personal care.

Designed as part of the wider Virtus seating range, the base integrates directly with the Virtus Seating System to provide a practical and supportive solution for educational environments. Its smooth manoeuvrability, substantial weight capacity and powered adjustment features help deliver long-term support for pupils, carers and school staff.

Key Features

- Powered height adjustment and tilt functionality
- Heavy-duty frame for everyday school use
- Supports learning, play and therapy activities
- Available with manual or powered adjustment
- Smooth manoeuvrability with lockable castors

▶ TECHNICAL SPECIFICATIONS

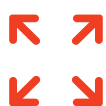

Base Size

Width: 65 cm

Depth: 83 cm



Weight Limit

User: Max. 90 kg

+ **Seat:** Max. 120 kg

Tilt-in Space

5° of Forward Tilt

25° of Backward Tilt

Height Adjust

37.5cm at lowest


62.5cm at highest



Adjustable Push Handle

Min. 70cm

Max. 90cm

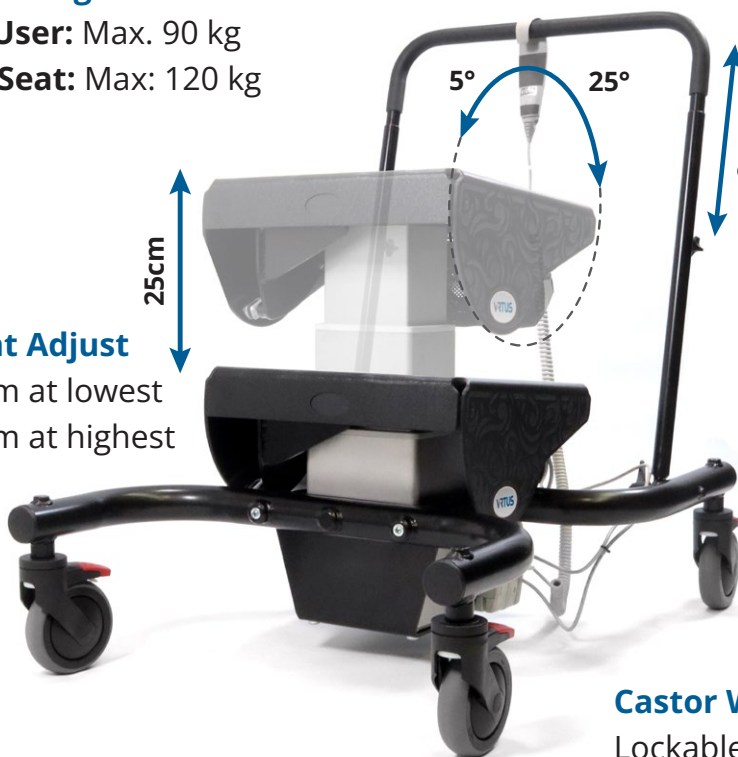

Castor Wheels

Lockable Ø125mm

Battery



- Up to 50 lift cycles per charge
- Can be operated whilst plugged in
- Integrated cable storage basket

POWERED ADJUSTMENT

The Virtus Powered High-Low Base provides smooth powered height adjustment with a 25cm lift range, supporting transfers, interaction and everyday positioning throughout the school day.

CLASSROOM PARTICIPATION

Helps position users at appropriate heights for classroom activities, group work, mealtimes and therapy sessions. Height and tilt adjustment allow seating positions to be adapted throughout the day, helping maximise comfort, participation and engagement.

POWERED TILT

Powered forward and rear tilt provides up to 5° forward tilt and 25° rearward tilt, supporting positioning, comfort and hoisted transfers whilst helping maintain user stability.

DESIGNED FOR SCHOOL LIFE

The robust steel frame combines smooth manoeuvrability, lockable castors and practical adjustment features to support the demands of busy educational environments.

MANUAL HIGH-LOW BASE

The Manual High-Low Base provides a simple and reliable solution for positioning children at appropriate heights for learning, play and everyday activities. Height adjustment allows the seating system to be adapted to different classroom environments whilst supporting interaction with peers and participation in group activities.

Its robust design, ease of use and lower acquisition cost make it a practical choice for many educational settings.



User + Seat Max. 60kg
-10° - 30° Tilt

Designed Around Educational Environments

Supporting Learning and Participation

- Suitable for classroom, dining and therapy activities
- Available with manual or powered High-Low Bases
- Supports interaction with peers and staff
- Designed to adapt alongside changing needs
- Extensive range of postural support options

Manufactured in the UK

- **ISO 13485 certified**
- **Durable components** designed for long-term use
- Responsive **technical support** and after-sales assistance



Choice of 7 machine washable covers



The Virtus seating system is designed to be robust and provide safe, effective postural support for an average lifetime of 8-years, under normal use conditions.

Every component on the Virtus can be replaced to allow for continued use and re-issue following servicing & maintenance.

The actual service life of individual products may vary depending on frequency of use, environment and maintenance. The service life can be maximised by adhering to the manufacturer's IFU.

Disclaimer: *The above applies if the product is under normal use conditions, free from unreasonable or deliberate misuse.*



Scan to
download
the Virtus
Instructions

BOOK A FREE ASSESSMENT

Our nationwide team of Product Specialists can provide **assessments, demonstrations** and **product support**, helping schools configure each Virtus Seating System to meet the individual needs of the user.





Rehabilitation Manufacturing Services

Thompson House, Unit 10, Styles Close,
Sittingbourne, Kent ME10 3BF

Tel: +44 (0) 1795 477280

Email: sales@rms-rehab.co.uk

www.rms-rehab.co.uk