





Juditta

An evolution in comfort.

Juditta is a **tilting and reclining seating solution** designed and built by Ormesa for teen and adult users who need a supportive seat for daily activities in offices and schools.

It is a seating system recommended for people with **different levels of disability or postural support**, even adapting to the most complex cases. It allows someone to achieve a **stable and comfortable position** and to prevent incorrect postural positioning.

The **viscoelastic chords** adapt to each individual user's anatomy and prevent rubbing and pressure sores. The particular structure allows **breathability, hygiene and comfort**.

PRODUCT
MARKED 

JUDITTA is a MEDICAL DEVICE that COMPLIES WITH THE SAFETY REQUIREMENTS defined with the REGULATION (EU) 2017/745

Hygienic and breathable

Anti-bacterial headrest

made of a special material with closed cells
highly resistant to external chemical agents.

Long-lasting visco-elastic material

The seat and backrest are covered with viscoelastic
fire-proof material. It is resistant, breathable and
anti-slip. It is made of a special and
exclusive Ormesa formula.



Ergonomic sitting

Ergonomic backrest

It holds the shapes of the user. The special structure of the backrest gives a proper alignment to the trunk and **it doesn't generate any pressure on the spine.**

It also offers a high level of stability without the need of side supports. It holds deformities and prevents retractions.

It relieves and redistributes the pressure and increases the sitting position tolerance

Backrest adjustable in height to adapt to the different dimensions of users.

Seat

The **special cover of the seat** and the slight tilt prevent the user from sliding forward.

The leg and backrest mechanical joints are close to the physiological ones. When you change the tilt of the leg rest and backrest, otherwise unavoidable slips are not created.

Multi-adjustable headrest

The structure and dimensions of the headrest provide stability to the user without obstructing the lateral movement of the head.

The headrest follows the shape of the backrest; thus assuring the proper alignment between head and trunk. With depth and lateral adjustment, the headrest can be offset from the alignment of the trunk and pelvis.

All the adjustments facilitate and sustain the visual orientation of the user.




**Lateral
Translation
Adjustment**



**Angle
Adjustment**



**Depth
Adjustment**

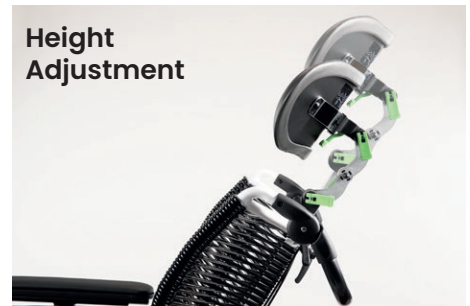

**Height
Adjustment**



Tilt-in space and reclining

The tilting and reclining features make it possible to get the **ideal support position**.

This allows a variation of the pressure points and prevents from generating pressure sores.

Tilting and recline are operated by two independent gas springs that absorb the user's weight, making adjustments easy and allowing postural settings to be changed multiple times throughout the day.



Comfortable

Distributed pressure discharge


Weight is well distributed between the seat and backrest, excluding pressure peaks on the ischial area.

The viscoelastic material of the seat and backrest follows the anatomy of the user adapting to it without generating any pressure and/or rubbing.

Juditta sitting shape and cover avoid pressures along the spine.




Without a seated person

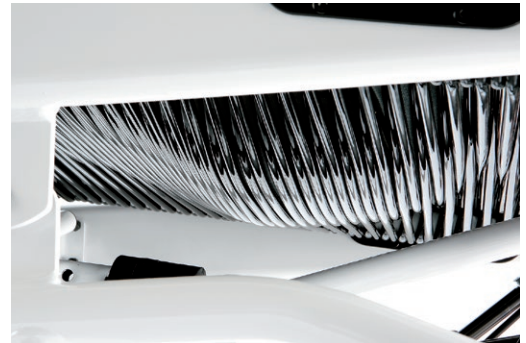

With a seated person

Clinical observation

Clinical observations carried out on Juditta show an optimal redistribution of the body weight. It entails the reduction of pressure peaks and prevents from the risk of generating pressure sores.

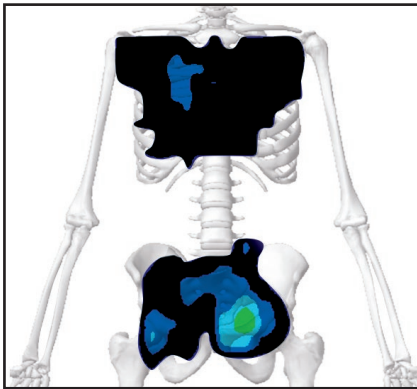


Juditta, our seating concept for those patients who need to maintain the lying position for a long time.

The special viscoelastic cords adapt themselves to the shapes and weight of the user, it increases backrest and seat contact area and generate a considerable reduction of the pressure peaks.

*Juditta: max. tilting, max. backrest recline, lifted footrests.
Patient (average adult) just positioned.*




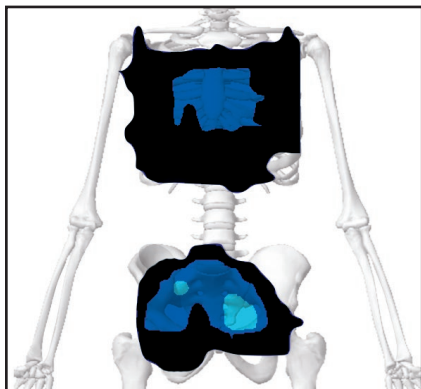

Extended pressure discharge

Even after one hour and 30 minutes from positioning the user the special viscoelastic cords continue to maintain their shaping feature and the pressure loading decreases.

*Juditta: max tilting, max backrest recline, lifted footrests.
User (average adult) positioned by 90 minutes.*

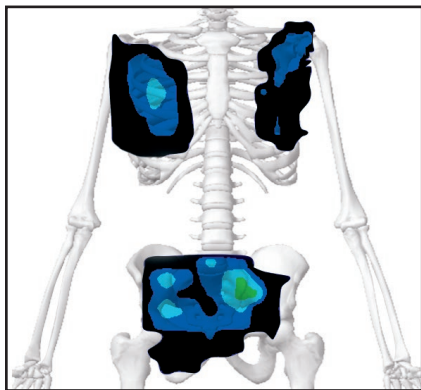

No pressure loading on the vertebral bodies

On the left we see a safety situation for the user with limited pressure discharge. When removing the upholstery there are no pressures along the spinal tract. This makes Juditta particularly suitable for people suffering from inflammation and skin changes in the area of the spine.

Juditta: results without upholstery.

No risk

Customisable

Adaptable to different user needs



Legrest

Removable, swing away with adjustable elevation without effort, activated with gas springs.

Calf rests adjustable in depth and free tilt to accommodate the movements of the leg.



Footrests

Adjustable in flexion-extension and height.

Customisable

Adaptable to different user needs



Armrests

Padded and easily adjustable in height.



Sides with scissor opening

A pivot joint behind the backrest allows the armrest to swing completely away, providing unobstructed access for transfers.



Seat depth adjustment

Four seat depth positions in 3 cm increments, while maintaining the same backrest recline angle at each setting.



Backrest recline

easily operated by levers.



Backrest

folding for transportation and storage.

Breathable, durable and long lasting upholstery

- **Fire-proof:** its fire-proof properties are inherent to the material and they do not come from any additional treatments.
- **Recyclable:** production cycle with low environmental impact.
- **Martindale method:** maximum resistance to rubbing: 100.000 cycles. UNI EN ISO 12947.
- **Oeko-Tex Certificate:** for all fabrics that are in direct contact with the skin, to guarantee quality and safety for health.
- **Long-lasting and highly resistant colours:** they are resistant to sweat, light, dry cleaning and washing machine (*see its Use and Instruction manual*).



CERTIFICATE
MED annex A. 1/3. 19-20-21



Accessories



Side-Temporo Support
Adjustable in Height & Depth



Side-Temporo-Fronto
Adjustable in Height & Depth



Removable Headrest Cover



45° Pelvic Belt



Padded Abduction Block
Adjustable in Height & Depth



Water-Proof Cover



Reclinable Tray
Tilt & Depth Adjustable



Tray Covering
Padded & Removable



IV Pole



Footplate Covering
Padded & Removable




Footrest lock



Padded Leg Panel
Alternative to the Calf Rests

Juditta B12


size and weight

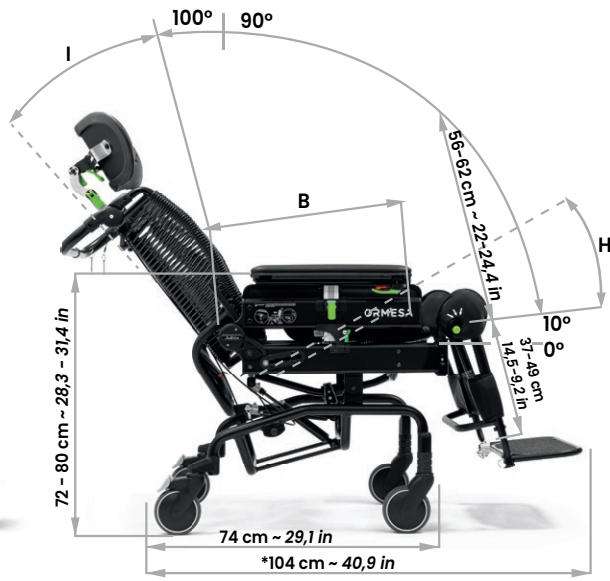

*footrest adjusted to the maximum length



	size 36	size 40	size 45	size 50
A: seat bottom width	36 cm ~ 14,1 in	40 cm ~ 15,7 in	45 cm ~ 17,7 in	50 cm ~ 19,6 in
B: seat bottom depth	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in	40-43-46-49 cm 15,7-16,9-18,1-19,2 in
C: backrest width	47 cm ~ 18,5 in	47 cm ~ 18,5 in	52 cm ~ 20,4 in	57 cm ~ 22,4 in
D: overall width	60 cm ~ 23,6 in	60 cm ~ 23,6 in	65 cm ~ 25,5 in	70 cm ~ 27,5 in
H: tilt-in space	from 10° to 35°	from 10° to 35°	from 10° to 35°	from 10° to 35°
I: backrest inclination	from 10° to 38°	from 10° to 38°	from 10° to 38°	from 10° to 38°
Total weight	37 kg ~ 81,5 lb	37 kg ~ 81,5 lb	39 kg ~ 85,9 lb	40 kg ~ 88,1 lb
Max capacity	120 kg ~ 264,5 lb	120 kg ~ 264,5 lb	130 kg ~ 286,6 lb	140 kg ~ 308,6 lb
User height	from 140 to 190 cm from 55,1 to 74,8 in	from 140 to 190 cm from 55,1 to 74,8 in	from 140 to 190 cm from 55,1 to 74,8 in	from 140 to 190 cm from 55,1 to 74,8 in

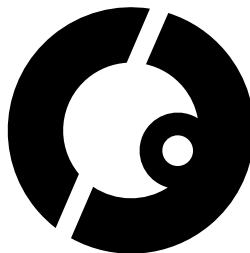

A special thanks to the families who participated with enthusiasm in our project. It is thanks to our customers and their stories that we find the energy and stimuli to evolve day by day.

COMPANY WITH QUALITY
SYSTEM CERTIFIED BY DNV
ISO 13485

Sold and distributed in the UK by



Rehabilitation Manufacturing Services



01795 477 280

www.rms-rehab.co.uk

sales@rms-rehab.co.uk

Thompson House, Unit 10 Styles Close,
Sittingbourne, Kent, ME10 3BF